



The *smarter source* for affordable specialty instruments

MASTEL DAVIS MD ONE-STEP MARKERS

PRE-OP MARKERS

- Designed for Toric IOL's, LRI and AK markings
- Rotating head of this plum-bob handle is laser engraved in 5° increments whereby interpolation allows setting within a degree or so
- Represents our best in applied precision and efficiency as the pre-marks themselves serve as the intraoperative toric alignment marks
- Eliminates secondary marking steps that add to potential translation errors
- **Fin model** - is simply three marks, two of which form the steep meridian and the third functions as a 90° cardinal reference point
- **Arc model** - features the steep meridian line as well as a partial capsulorhexis OZM
- **Ahmed model** - adds in arcs at a 9.5 mm OZM and marks at 90°, 75°, 45°, and 30°



Fin
MRMT-OF



Arc
MRMT-OA




Ahmed
MRMT-OI



Ambler Surgical is an authorized distributor for the Mastel Precision line of ophthalmic instruments and diamond scalpels. Mastel Precision has been involved exclusively in the design and production of ophthalmic instruments since 1981. Mastel Precision offers unparalleled design, quality and precision-crafted surgical instruments and diamond scalpels for astigmatism management and corneal procedures that are considered unmatched in the ophthalmic industry.

To view pricing and our complete line of Mastel Precision Instruments, visit us online.

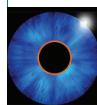

OPHTHALMOLOGY

amblersurgical.com/mastel

1-888-407-0006